


# ● THERMOGARD HQ™ Setup Guide


## 1 Insert Coil

- A. Turn on console
- B. Open cap
- C. Place coil in well
- D. Close cap, passing tubes through slit in cap

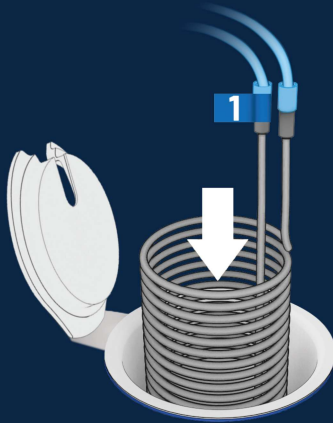

## 2 Place Air Trap

Place in stand with **orange** end down

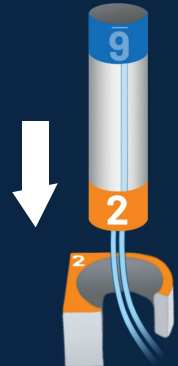

## 3 Load Pump Tube

- A. Open roller pump lid
- B. Fully insert tubing disc **into slot 3**
- C. Firmly rotate counter-clockwise while feeding tube
- D. Fully seat tube **into groove 3**



## 4 Spike Saline Bag

Ensure spike fully extends into saline bag



## 5 Prime Tubing

- A. Close roller pump lid
- B. Press and release button
- C. Pump will run for 2 min. Confirm that pinwheel turns and bubbles are cleared.

⚠ If bubbles continue to circulate after 2 min, recheck steps 2 – 5



## 6 Insert Air Trap

- A. Insert into chamber with **blue** end down
- B. Close console top cover

⚠ Ensure tubing is not pinched by top cover. See back side of page for proper placement.



# THERMOGARD HQ™ Setup Guide ●

Thick tubing  
in roller  
pump

Saline bag tubing  
not pinched by  
top cover

Saline bag tubing  
into slot

Extension tubing  
in notches


T1: Primary patient  
temperature probe

T2: Secondary patient  
temperature probe  
(optional)


**ZOLL**

ZOLL Circulation, Inc.  
2000 Ringwood Avenue  
San Jose, CA 95131 USA  
+1-877-225-7487  
[www.zoll.com](http://www.zoll.com)

Catheter connection for Treatment

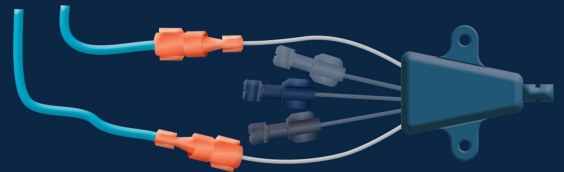

Catheter connection for Transporting

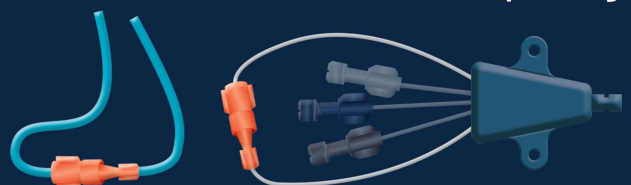

# ● THERMOGARD HQ™ Display Guide



# THERMOGARD HQ™ Display Guide ●

| 1               | 2            | 3                  | 4                   |
|-----------------|--------------|--------------------|---------------------|
| WiFi Connection | WiFi Battery | Digital EMR Output | Patient Temp Output |
| ● No Connection | ● Power Off  | ● Disabled         | ● Disabled          |
| ● Connected     | ● Charged    | ● Enabled          | ● Enabled           |
| ● Attn Required | ● Replace    | ● Attn Required    | ● Attn Required     |

5

## USB Data Download:

1. Plug in USB drive
2. Wait for solid green light
3. Remove USB drive

- No Connection
- Copying
- Finished
- Attn Required

