


Zenix® Viewer

Hospital Technical Specifications

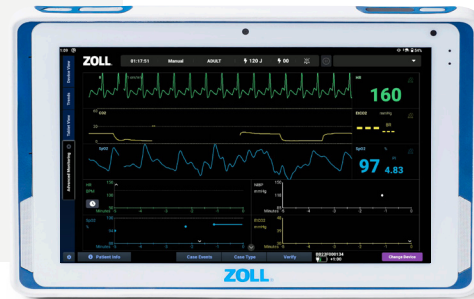

Zenix transforms the way your code leader or CPR coach provides guidance during a code with the first-of-its-kind Zenix Viewer. Without disrupting care, anyone holding the Viewer can be anywhere in the room—out of the way and fully informed with the exact data they need. Instantly customizable views include CPR with compression and ventilation data, vital sign trends, and playback and review.

The Viewer connects with RescueNet® CodeWriter for code documentation and RescueNet CaseReview for debriefing. The Viewer even helps with training, allowing educators to run mock codes from a distance.

System

Processor: Qualcomm SDA660 / SDM660 (up to 2.2GHz)

Memory: 4 GB LPDDR4x

Operative System: Android 13

Display

Type: 10.1" FHD LCD

Resolution: WUXGA 1920x1200

Brightness: 380 cd/m²

Touch Type: 10-point multi-touch P-CAP touch control

Glass Type:

Corning® Gorilla® Glass 3

Touch Mode: Normal, Glove

Communication

WLAN: IEEE 802.11 a/b/g/n/ac/k/v/r

Bluetooth: Bluetooth® v5.0 + BLE

NFC: 13.56 MHz, compatible with ISO15693, ISO14443A, ISO14443B, FeliCa

Standard: 1 x USB3.1 Type C supporting data transmission only*, 1x USB 2.0 Type A, 1x Combo audio jack, 1 x 19 VDC-in jack, 1 x SIM card reader, 1 x micro SD card

Extended: 1 x AIM extension 14-pin pogo connector: 1 x USB 2.0, 1 x 2-wire UART, 5V/1A power

Barcode Scanner (optional)

Integrated 1D/2D barcode scanner

Sensors: Ambient light sensor, e-Compass, G-sensor, gyroscope

Security

TPM 2.0

Warranty

2 years

ZOLL®

Power System

Input Voltage: DC jack: 19V/3.43A

Battery (Removable)

Smart battery with capacity LED indicators; 10.8V, 2400 mAh, 26 WHr

Operating Time:

Approximately 6 hr (depending on usage scenario)

Charging Time: < 3 hours

System Indicator Full:

Full: green (100%); Charging: blue; Low: red ($\leq 10\%$)

Environment

Operating Temperature:

0 ~ 40°C (32°F - 104°F)

Charging: 0 - 45 °C (32 - 113 °F)

Storage Temperature:

-20 ~ 60°C (-4 ~ 140°F)

IP Rating: IP65

Shock: Operating: 10 ~ 20 G, 11 ms; Non-operating: 50 G, 11 ms

Vibration:

Operating:

(5 ~ 500 Hz) 0.5 ~ 1 G;

Non-operating:

(5 ~ 500 Hz) 0.04 PSD

Drop-tested: Up to 120 cm (4 ft) onto plywood over concrete

Humidity Operating: =10 ~ 90%

@ 40 °C non-condensing

Humidity Storage: 10 ~ 95%

non-condensing

Physical Characteristics

Dimensions: (W x D x H) 280 x 18.3 x 180.1 mm (11 x 0.7 x 7.1 in)

Weight: 980 g (2.2 lb)

Color: White/blue

Control Buttons:

1 x power button, 2 x volume key (up/down)

Certifications: IEC/UL/EN 60601-1, FCC/IC, CE, NCC/BSMI, TELEC/VCCI, RCM, CCC/SRRC

ZOLL Medical Corporation

269 Mill Road | Chelmsford, MA 01824 | 978-421-9655 | 800-804-4356 | zoll.com

Specifications subject to change without notice.
Copyright © 2025 ZOLL Medical Corporation. All rights reserved. Zenix and ZOLL are trademarks or registered trademarks of ZOLL Medical Corporation in the United States and/or other countries. All other trademarks are the property of their respective owners.
Printed in U.S.A. MCN HP 2510 0763

The ZOLL logo consists of the word "ZOLL" in a bold, white, sans-serif font, centered within a solid blue rectangular background.